

# 5mL Magnetic Tube Dock (MTD-5)

Model # A4103

**WAVESENSE**

## **Intended Use:**

For In Vitro Diagnostic Use.

This product is intended for the recovery/enrichment of magnetic particles from biological fluids and tissue culture.

## **Description:**


The MTD-5 is a reusable magnetic dock containing a single powerful neodymium permanent magnet suitable for isolating submicron and micron paramagnetic particles. It can be used whenever it is desirable to debulk specimens laden with debris, bacteria, mucous, and dense background cell populations typically found in cell pellets obtained from urine, blood, bone marrow, and tissue homogenates.



## **Supplied As:**


Catalog #	Contains
A4103-1	1 Each
A4103-5	5 Pack

## **Featured Components:**

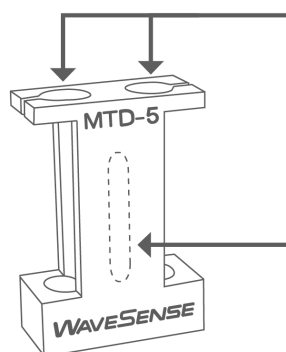

Channeled design firmly aligns 12 x 75 mm culture tubes with the magnetic field.

Embedded magnets attract magnetically labeled antibodies to the side or base of the tube.

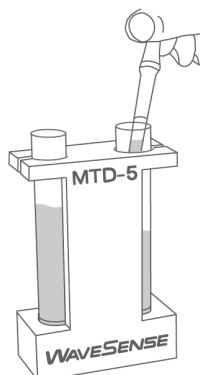

Illustrated with 12 x 75 mm culture tubes and pipette (not provided).

## **Storage:**

Store dry at room temperature.

## **Compatible Accessories (not included):**

- Standard 12 x 75 mm culture tubes (available from multiple suppliers)

## **Material Safety Data:**

Contains plastic components and magnets. In the event of breakage, wear gloves and eye protection for clean up as some pieces may be sharp and could cause injury. Keep away from electronics or cards with magnetic strips.

Product Specification Sheet



15339 Barranca Pkwy  
Irvine, CA 92618 USA  
www.WaveSense.net

Toll Free: 800.807.7760  
Phone: 949.341.1980  
Fax: 949.341.1982  
Contact@WaveSense.net

ML111128-2  
© 02/01/12

Page 1 of 1